

1. AT-22 tablets: stability review

Aldreversa Therapeutics | WA/2024/211 | 25 Jul 2025

2. Approved numerical evidence

Condition	Timepoint (months)	Assay (% LC)	Total impurities (%)	Result
25 °C/60% RH	0	99.7	0.17	Complies
25 °C/60% RH	1	100.1	0.14	Complies
25 °C/60% RH	3	99.7	0.22	Complies
25 °C/60% RH	6	99.8	0.22	Complies
25 °C/60% RH	9	99.1	0.31	Complies
25 °C/60% RH	12	99.7	0.32	Complies
30 °C/65% RH	0	99.7	0.17	Complies
30 °C/65% RH	1	100.5	0.22	Complies
30 °C/65% RH	3	99.5	0.24	Complies
30 °C/65% RH	6	99.4	0.31	Complies
30 °C/65% RH	9	98.8	0.35	Complies
30 °C/65% RH	12	98.9	0.44	Complies
40 °C/75% RH	0	99.7	0.17	Complies
40 °C/75% RH	1	99.4	0.22	Complies
40 °C/75% RH	3	99.4	0.42	Complies
40 °C/75% RH	6	97.8	0.67	Complies


3. Protocol identity and sample custody

- Client: Aldrevessa Therapeutics
- Protocol: WA/2024/211 | LIMS: WST-24-006564
- Opened 22 Jul 2024; completed 22 Jul 2025
- Dr L. Wu: stability lot and sample-custody reconciliation.
- Storage condition and pull interval are carried with the stability lot on the batch card.
- Sequence sheets retain the pack or presentation identity; results must not be pooled across storage conditions.
- Tested interval: through 12 months.

4. HPLC method qualification

Parameter	Result	Acceptance
Linearity (r^2)	0.9991	≥ 0.999
Range (% nominal)	50-150	80-120 minimum
Accuracy 80% (% recovery)	99.4	98.0-102.0
Accuracy 100% (% recovery)	99.4	98.0-102.0
Accuracy 120% (% recovery)	99.6	98.0-102.0
Repeatability (%RSD, n=6)	0.41	≤ 2.0
Intermediate precision (%RSD)	0.72	≤ 2.0
LOQ (% of nominal)	0.05	≤ 0.05

5. Impurity trend by storage interval

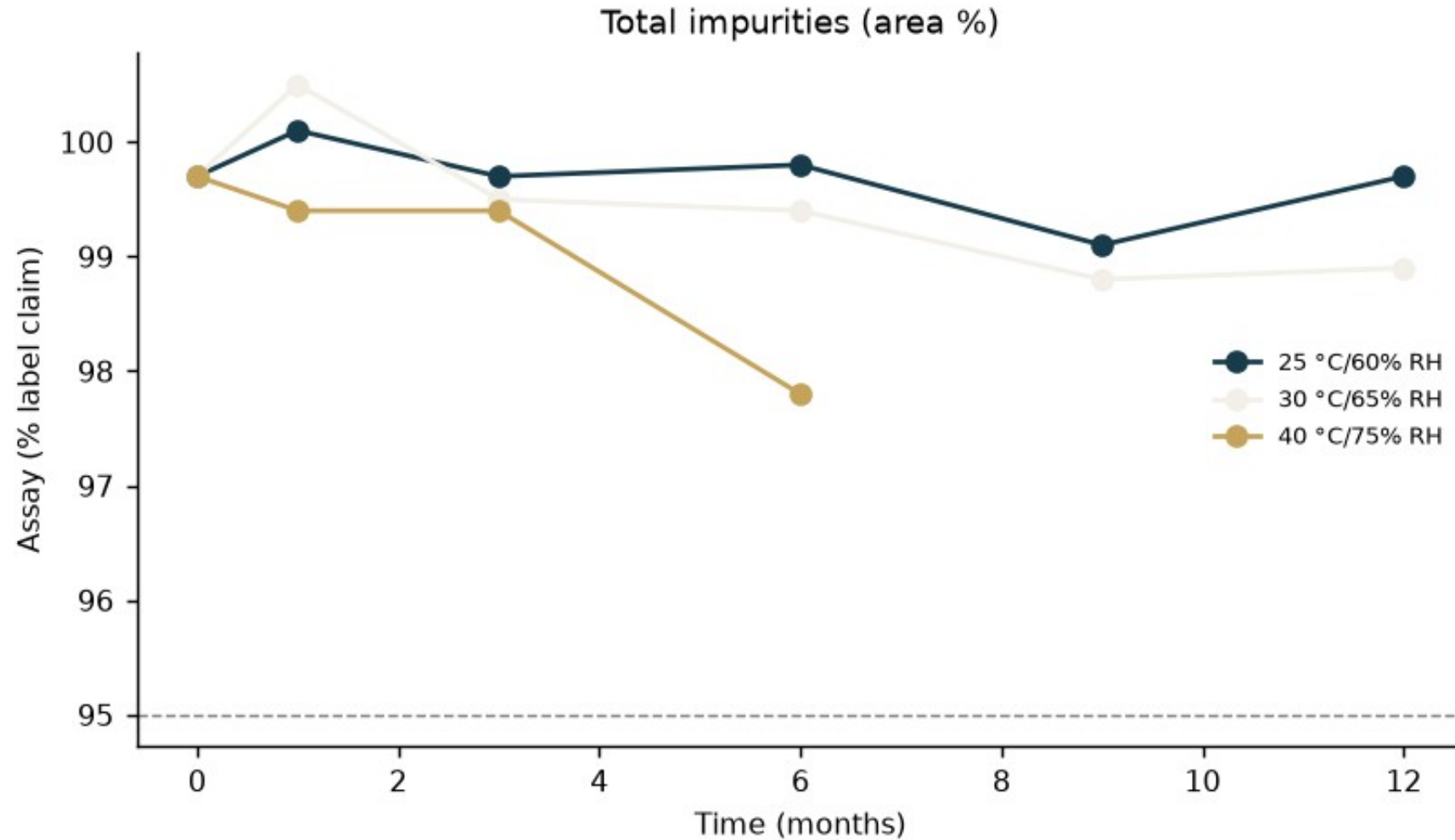

6. System suitability and deviation disposition

- Saira Mahmood: analytical method and system suitability review.
- Dr L. Wu: reconciliation of preparation, storage and sequence identity.
- Dr A. Penrose: approved interpretation and client release.
- Approved sequence sheets are the authority for suitability acceptance and any repeat sequence.
- Deviation notes remain linked to the affected batch card; this deck does not replace their disposition.
- No additional chamber excursion, reinjection or acceptance limit is inferred from the trend plot.

7. Interpretation for transfer

- The uncoated comparator showed a humidity-driven impurity increase.
- Aldrevessa selected the coated tablet.
- Carry the coated tablet into transfer; retain the uncoated comparator results as the formulation decision record.
- Conclusion applies to tested storage conditions and the reported 12-month interval.

8. Release gate and client decision

- Aldrevessa selected the coated tablet.
- Release basis: AT22_PULLS_AND_SUITABILITY_RESULTS_Rev03.xlsx; approved 22 Jul 2025.
- Carry the coated tablet into transfer; retain the uncoated comparator results as the formulation decision record.
- Archive the approved readout with the sequence sheets, batch cards and linked deviation notes.